

# Zirkonzahn<sup>®</sup>

*Human Zirconium Technology*




**ALL ZIRCONIA**

## *Pre-sinter Colouring*

Milling procedure complete, the bridge frame is cut from the block. Labial contours are refined, but a support base for sintering remains.

The frame is coloured in various shades with Zirkonzahn Colour Liquids. Extra care is taken in defining the border between dentin and tissue areas. After colouring the frame is dried under an infra-red heat lamp for 60 minutes then placed in the sinter furnace in vertical position.

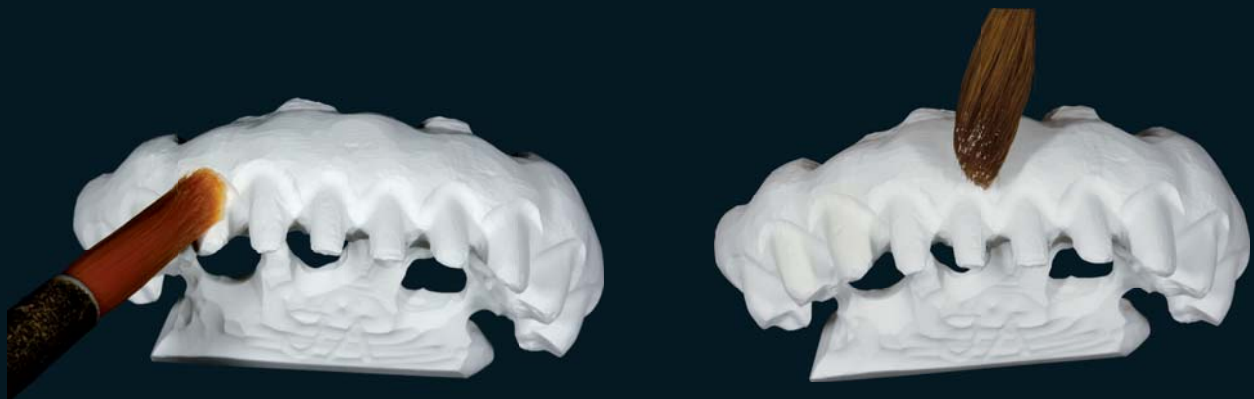

## *Result after Sintering*

Master Technician Juergen Auffarth, Salzburg – Austria



# Zirkonzahn®



## HUMAN ZIRCONIUM TECHNOLOGY

ZIRKONZAHN USA Inc. - Tel 1 800 989 8931 - Tel 678 441 9419 - [infousa@zirkonzahn.com](mailto:infousa@zirkonzahn.com)

ZIRKONZAHN World Wide - Tel +39 0474 066 680 - [info@zirkonzahn.com](mailto:info@zirkonzahn.com) - [www.zirkonzahn.com](http://www.zirkonzahn.com)



*Milled with the Zirkonzahn Milling System*