

# Zirkonzahn<sup>®</sup>

*Human Zirconium Technology*



**100% PRETTAU ZIRCONIA**

## *A lower patient-removable full denture – guaranteed to last*

The idea was to produce a bar-born removable lower restoration made entirely from ceramic material. Two zirconia bars were milled; both fitted with slide-attachments, one at each end. Through Teflon-sleeve matrixes over the attachments sufficient friction is achieved to hold the denture in place and give the desired stability. Such longevity is guaranteed.


The bars are milled by reducing a full-anatomical-contour mock-up frame; in this way optimal bar position within the restoration is ensured. The entire prosthesis is made from 100% full zirconia. Only the pink tissue was layered with veneer porcelain. The result is a Prettau Bridge, worn over bars and entirely patient-removable.



## *The Pretttau Bridge*

Georg Walcher, Zirkonzahn Italy

Dr. Nico Bühler – Clinic for Reconstructive Dentistry – University Basel/Switzerland

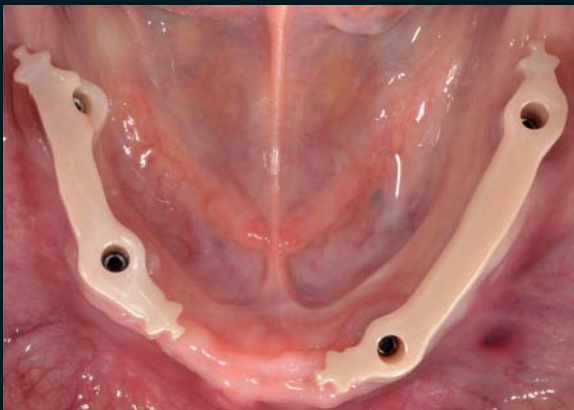

# Zirkonzahn®



## HUMAN ZIRCONIUM TECHNOLOGY

ZIRKONZAHN USA Inc. - Tel 1 800 989 8931 - Tel 678 441 9419 - [infousa@zirkonzahn.com](mailto:infousa@zirkonzahn.com)

ZIRKONZAHN World Wide - Tel +39 0474 066 680 - [info@zirkonzahn.com](mailto:info@zirkonzahn.com) - [www.zirkonzahn.com](http://www.zirkonzahn.com)



WEAA5409=

*Milled with the Zirkonzahn Milling System*